
Kegel Master


Gen 2

Your Personal Coach Kegel Exercise

EC REP Interslash BV
Tomatenmarkt 1, 1681 PH Zwaagdijk, NL

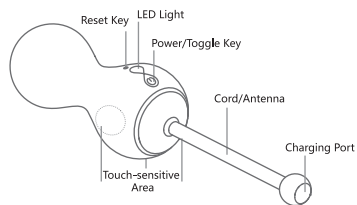Package List

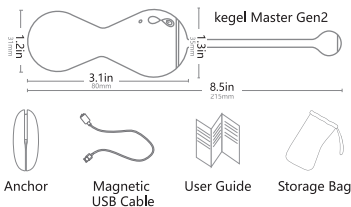Specification

Material: Body-safe silicon /ABS | Weight: 58g | Battery: 200mAh
 Charging voltage: 5V | Maximum Charging Current: 500mA
 Charging time: 60 minutes | Working time: 180 minutes
 Waterproof: IPX6

Smartphone Workout Mode

Step1: APP Download

1. Search for "Magic Kegel" in IOS App Store or Google Play
2. Scan below QR code to download directly.
3. Download from www.magicmotion.fun


Step2: Run "Magic Kegel" APP


Turn on your Kegel Master Gen 2 device

Follow App guide to connect

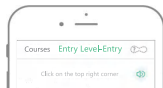It's connected successfully when it shows "connected"

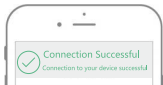Step3: Enjoy Kegel Training

Find a comfortable position for inserting. Many women find lying on the back or standing as most comfortable. Kegel Master Gen 2 LED facing upward, please check below picture. For pushing forward easily, apply a bit water-based lubricant on the top of device or on the condom.

Tips: Please leave the cord of Kegel Master Gen 2 outside the body, at least 4 cm length, to keep good signal without barrier.

Step4: Power Off

The device will turn off after pressing Power Key for 3 seconds, the LED light will turn off to indicate turn-off status.


Tips: Kegel Master Gen 2 will turn off automatically when it disconnected with App in 5 minutes.

Independent Mode (Without SmartPhone)


- Squeeze for vibration

Turn on Kegel Master Gen 2 and squeeze pelvic floor muscles to control the strength of vibration of device.

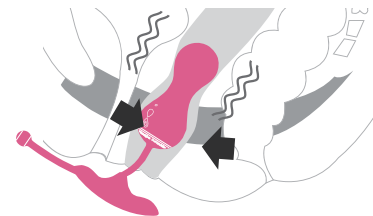
Charging Battery

Step1 Open the cap of the magnetic charging port for charging.

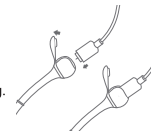Step2 Connect USB port to smartphone power adapter, laptop, Computer or power charging bank.

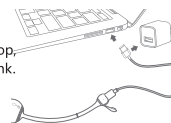Maintenance

Step3 LED light on to show charging normally. The LED light will be off when the charging is finished.

Tips: For first use, please charge the battery for 1 hour after the power is fully used.

Cleaning Magic kegel Master Gen2

Before Usage


Please clean Kegel Master Gen 2 with cleanser or warm water and anti-bacterial soap before first use.

After Usage

Please clean Kegel Master Gen 2 with cleanser or warm water and anti-bacterial soap and dry it.

Storage

Keep Master in a cool and ventilated place. Avoid direct sunlight and keep it away from children.

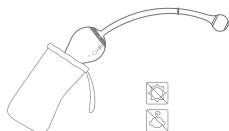Warranty Terms & Conditions

Magic motion warrants this product for a period of ONE (1) YEAR, after the date of original purchase, against defects due to faulty workmanship or materials. If you discover a defect and notify magic motion during the warranty period, magic motion will, at its discretion, replace the product free of charge.

The warranty covers working parts that affect the function of the product. It does NOT cover cosmetic deterioration caused by fair wear and tear or damage caused by accident, misuse or neglect. Any attempt to open or take apart the product (or its accessories) will void the warranty.

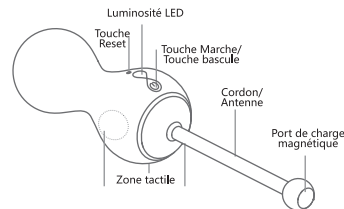
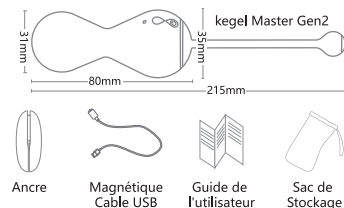
Claims under warranty must be supported by reasonable evidence that the date of the claim is within the warranty period. To validate your warranty, please keep your original purchase receipt together with these warranty conditions for the duration of the warranty period.

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